


Single Use Mechanical Lithotripter

LithoCrush V™

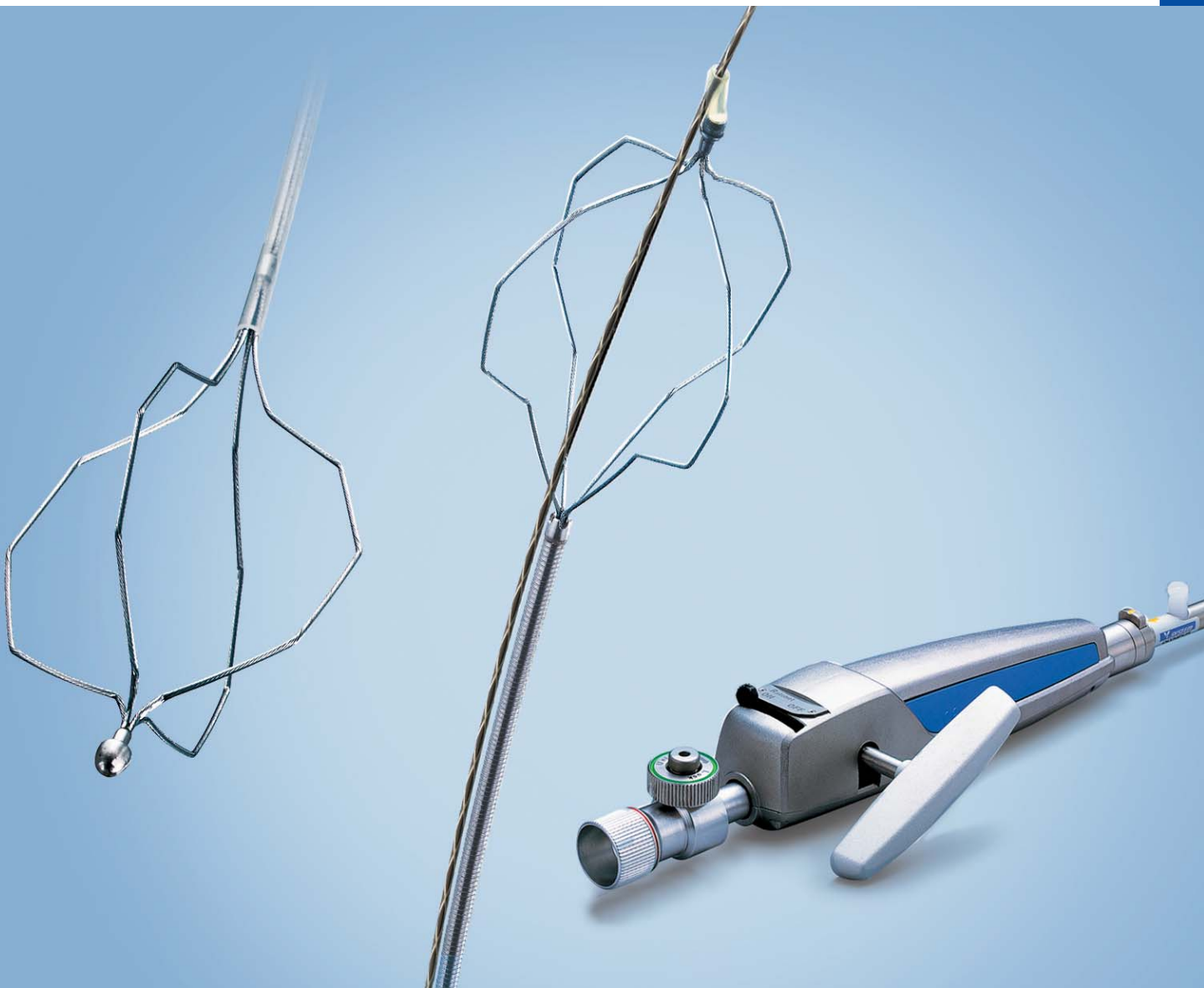


The easy-to-use, rotatable LithoCrushV™ features distal wireguided and bullet tip design options to fit the needs of virtually every procedure

LithoCrushV™


Bullet tip design provides powerful crushing performance

With a dual tube and coil sheath construction, the bullet tip design combines outstanding insertion ability with remarkable torque performance.

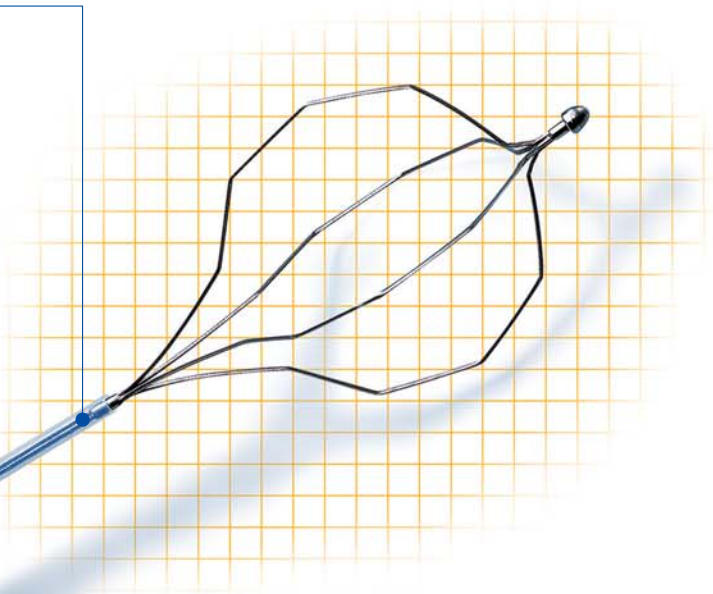

Extra strong basket available for enhanced crushing power

The hard-wire type basket with 19 strands versus the standard 12 strands has greater crushing power and is ideal for repeated crushing.

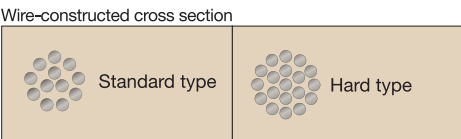



Convenient injection port

An injection port allows for injection of contrast medium with a syringe.



Distal wireguided design improves insertion capability

Distal wireguided design makes it easier to exchange catheters and insert the basket into the papilla. It also facilitates selective cannulation of the intrahepatic bile ducts or other target sites.

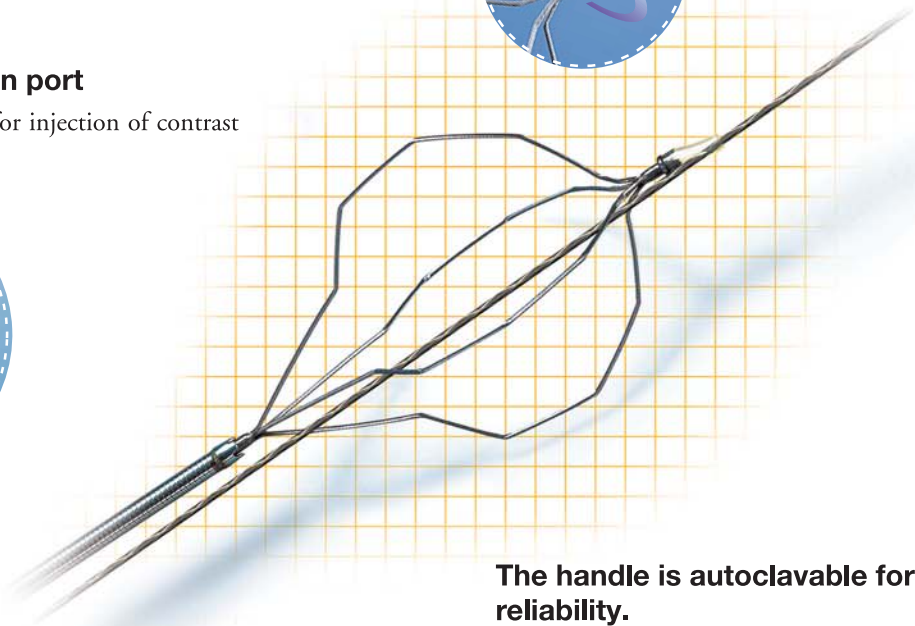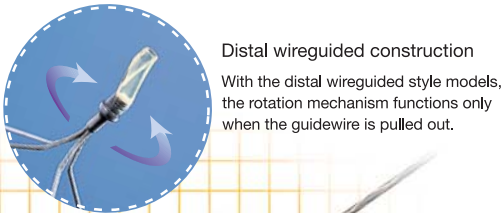


The handle is autoclavable for ensured reliability.

New handle design increases torque performance and simplifies setup

Advanced rotation mechanism improves retrieval performance

All models now feature a basket rotation mechanism that makes retrieval of small stones and crushed stones easier and more efficient.

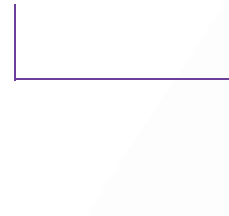


Ratchet prevents loss of traction on the stone

The ratchet prevents counter-rotation that can occur when a large, hard calculus is crushed.

Easy setup

The basket clicks into place effortlessly.

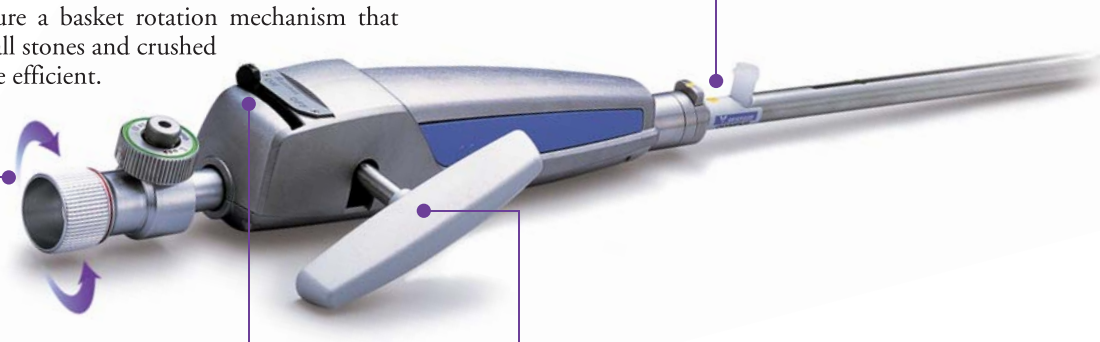

Increased torque ability

The contoured handle and T-Lever are easy to grip making it easier to apply pressure and control the basket when a stone is crushed.

Single Use Mechanical Lithotripter



Specifications

Model	Bullet tip/Wireguided type	Opening width	Hard/Standard basket-wire type	Minimum channel size	Compatible guidewire
BML-V242QR-30	Bullet tip type	30 mm	Hard	4.2 mm	-
BML-V237QR-30	Bullet tip type	30 mm	Hard	3.7 mm	-
BML-V232QR-30	Bullet tip type	30 mm	Standard	3.2 mm	-
BML-V232QR-26	Bullet tip type	26 mm	Standard	3.2 mm	-
BML-V437QR-30	Distal wireguided type	30 mm	Standard	3.7 mm	0.89 mm (0.035")
BML-V442QR-30	Distal wireguided type	30 mm	Hard	4.2 mm	0.89 mm (0.035")

LithoCrushV handle

: Refer to instructions.

Model	Feature
MAJ-441	Rotation mechanism, ratchet prevents counter-rotation, autoclavable



BML-110A-1

In the unlikely event of basket impaction, the BML-110A-1 Emergency Lithotripter Handle is to be used.

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.



OLYMPUS MEDICAL SYSTEMS CORP.
Shinjuku Monolith, 2-3-1 Nishi-Shinjuku, Shinjuku-ku, Tokyo 163-0914, Japan

For a complete listing of
sales and distribution locations visit:
www.olympus.com