


Your Vision, Our Future

EVIS EXERA III

EVIS EXERA III GASTROINTESTINAL VIDEOSCOPE

GIF-1TH190

Therapeutic endoscope with superb HDTV image quality and a 3.7 mm instrument channel

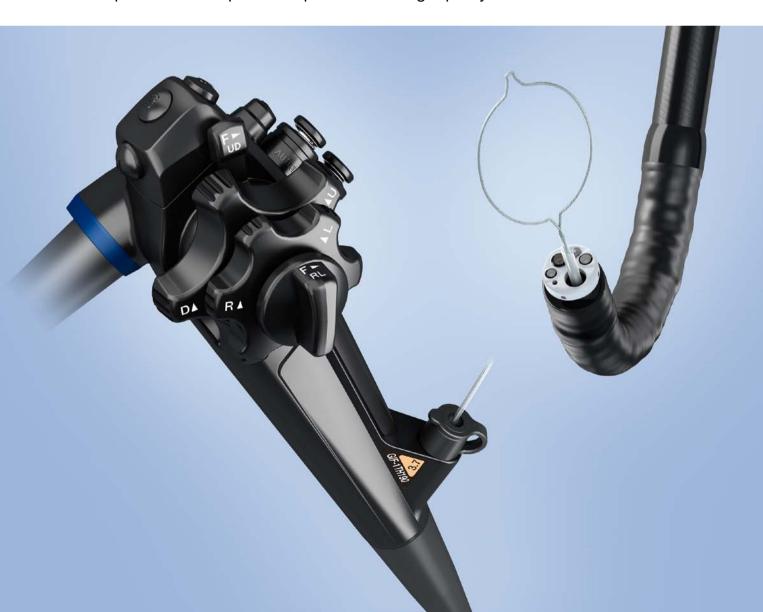

OLYMPUS GIF-1TH190



Main Features

HDTV Image Quality

HDTV image quality delivers high-definition observation capabilities, even in a new, slimmer scope design.



NBI (Narrow Band Imaging)

NBI in the GIF-1TH190 provides twice the viewable distance of conventional EVIS EXERA II 180 series scope*1 and is significantly brighter. HD NBI provides contrast, which may aid in the interpretation of mucosal morphology, vascular patterns, and blood vessel appearance.



Slim Design and Therapeutic Capability

The 10.0 mm distal end diameter is 0.9 mm smaller than the predecessor*2 – while still maintaining the flexibility of a 3.7 mm instrument channel diameter.

Close Focus

The GIF-1TH190 enables close observation of the mucosal surface (2 mm) without the need for electronic magnification.

Water Jet

A forward water jet helps to keep the mucosa clear during observation and treatment.





Waterproof One-touch Connector

A new connector design minimizes the effort required for setup prior to and in between cases. In addition, it is fully submersible and eliminates the need for a water-resistant cap and the associated risk of an expensive



3.7 mm

repair due to accidental immersion.



Specifications

Optical System	Field of view	140°	
	Direction of view	Forward viewing	
	Depth of field	2-100 mm	
Insertion Section	Distal end outer diame	neter 10.0 mm	
	Distal end enlarged		
	Objective lens	Up	Air/Water Nozzle
		XOX	Light-guide Lens
	Right	(00) [.eft
	Light-guide Lens Auxiliary Water Channel	Ins	trument Channel Outlet
	Down		
	Insertion tube outer diameter 10.9 mm		
	Insertion tube outer diameter 10.9 mm		
	Working length	1030	mm

^{*1} EVIS EXERA II Gastrointestinal Videoscope OLYMPUS GIF TYPE Q180
*2 EVIS EXERA Gastrointestinal Videoscope OLYMPUS GIF TYPE 1TQ160

Minimum visible distance 3.0 mm from the distal end Instrument Channel Direction from which endotherapy accessories enter and exit the endoscopic image Up 210° Down 90° Bending Section Angulation range Right 100° Left 100° Total Length 1350 mm Compatible EVIS Video System Center OLYMPUS CV-190 **EXERA System** Xenon Light Source OLYMPUS CLV-190

Channel inner diameter

Image courtesy of Robert H. Hawes, MD and Muhammad Hasan, MD (Florida Hospital Orlando, US) Water Jet image: Image courtesy of Horst Neuhaus, MD (Evangelisches Krankenhaus Düsseldorf, Germany)

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

